

exa-PACS · exa-RIS

Feature Summary

Merging Patient Charts


© 2023 Konica Minolta Healthcare Americas, Inc.

Merge patient charts


1. On the navigation (“burger”)  menu, select **PATIENT**.



2. Search for both patient charts by typing patient demographics as criteria.

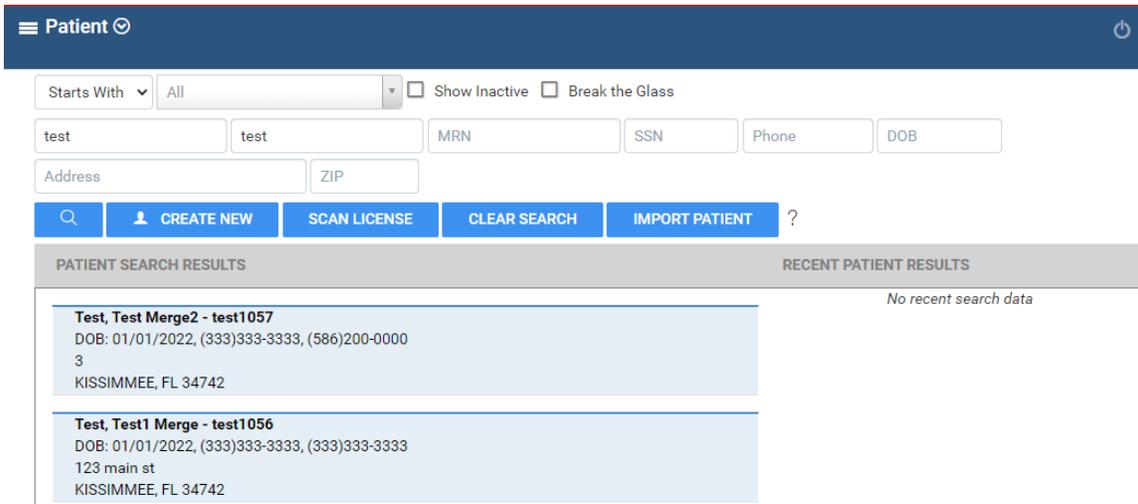

3. On the right side of the row with the *incorrect* patient chart, hover and select the gear button, and then select **Merge**.

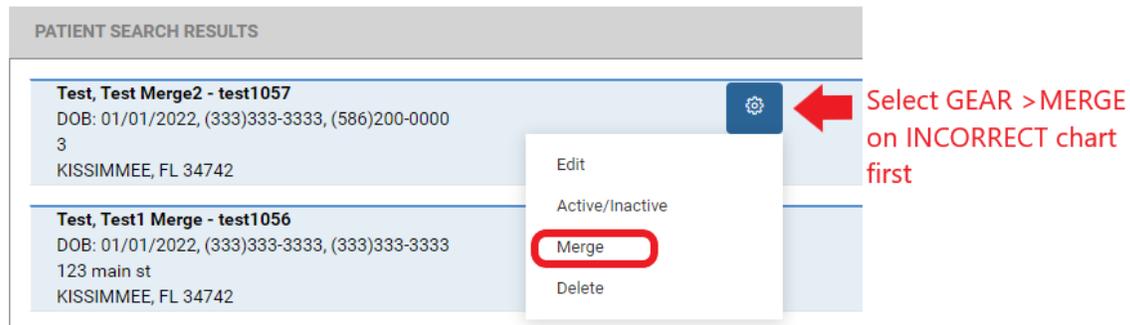

Result: A message prompts you for the *correct* patient chart.

Test, Test Merge2 Selected. Select second patient to merge into.

CANCEL

- On the right side of the row with the *correct* patient chart, hover and select the gear button, and then select **Merge**.
- In the **Patient Merge** screen, confirm the contents of the merge and then select **MERGE**.

Patient Merge
✕

Test, Test Merge2 (test1057)

01/01/2022 , M
3
KISSIMMEE, FL, 34742

➔

Test, Test1 Merge (test1056)

01/01/2022 , M
123 main st
KISSIMMEE, FL, 34742

MERGE

Description	Facility	Date
CT LEFT HUMERUS	Konica Minolta Healthcare Americas	Pre-Order
BONE DENSITOMET	Hodorville	Pre-Order
Chest X-Ray	Clemons OLD Facility	Pre-Order
COMPUTER-AIDED C	Hodorville	01/20/2023
DIAGNOSTIC MAMM	Hodorville	01/16/2023
DOP ECHO FTL PLST	DorothyLand	07/27/2022

Description	Facility	Date
Testing Different Mo	Hodorville	Pre-Order
Testing Different Mo	Pineapple Under The Sea	12/21/2022
Testing Different Mo	Pineapple Under The Sea	12/20/2022
Testing Different Mo	Pineapple Under The Sea	12/20/2022
Medial Branch Block	Clemons OLD Facility	12/13/2022
CT Thoracic w/o	Dave's Hospital	12/12/2022
XR CERVICAL SPINE	Dave's Hospital	12/12/2022
CT Cervical w/o	Dave's Hospital	12/12/2022
CT ABD C-/C+	Dave's Hospital	10/19/2022
CT ABD C-MATRL	DorothyLand	08/02/2022
DIAGNOSTIC MAMM	DorothyLand	07/29/2022
BONE MARROW BLC	DorothyLand	07/28/2022

Document type	Document Level	Created Date
No Records Found		

Document type	Document Level	Created Date
No Records Found		

Coverage Level	Insurance Name	Subscriber Info		
		Relation	Name	Address
Primary	AARP MEDI	Self	Test, Test	, , , ,

Coverage Level	Insurance Name	Subscriber Info		
		Relation	Name	Address
Primary	CIGNA	Self	Test, Test1	123 main st

Provider Name	Study Description
No Records Found	

Provider Name	Study Description
No Records Found	