



KONICA MINOLTA

SONIMAGE P3

Personal, Portable, Point-of-Care Ultrasound

Konica Minolta Medical Imaging introduces the new SONIMAGE P3, a personal, portable, point-of-care ultrasound system that revolutionizes Veterinarian practices and how you see your patients.



Giving Shape to Ideas

Adding Ultrasound to Your Veterinary

PERSONAL

The SONIMAGE P3 gives veterinary practices the ability to easily make use of ultrasound in their offices. Because of its smart-phone design, clinicians can focus on patients and scanning, not on cumbersome buttons and instructions.

PORTABLE

Being small and lightweight, the SONIMAGE P3 gives users the ability to quickly image at the point-of-care. Use as a stand-alone system or simply plug the hand held display unit or ultrasound probe directly into your Windows based PC, laptop or tablet running SONIMAGE VIEW software standard with the SONIMAGE system. It's flexible solution for occasions when you need a larger display.

POINT-OF-CARE

The SONIMAGE P3 enable veterinary offices the ability to offer ultrasound services, from routine checkups and procedures to emergency situations that simply require a look below the surface. Having ultrasound capability strengthens clinical confidence, aids in quick diagnosis, and enhances the patient and healthcare professional experience.

BROAD RANGE OF APPLICATIONS

- Check for free fluid – FAST exams
- Pregnancy checks
- Bladder assessments, volume calculations and assist with cystocentesis procedures
- Internal organ scans for gross abnormalities
- Trauma situations
- Lump analysis – cyst or solid
- Pleural effusions and pneumothorax
- Basic cardiac imaging



Dog Chest

“Our workflow has expanded by adding ultrasound to the daily practice. I do not need to refer cases to specialists that need basic sonograms as part of their work-ups, which benefits both our practice and the patient.”

Dr. Steven Smith, DVM

Practice Just Got Easier.

ANYTIME AND ANYWHERE YOU NEED IT

Whether you need to add a second machine to your practice or are looking to offer ultrasound to your services, the SONIMAGE P3 has a role to play. From sleep mode, it powers up instantly and is ready to go. Rather than having a dedicated ultrasound exam room, this device easily moves to the patient. Immediate ultrasound helps speed diagnosis, reduce patient wait-times, and improves workflow by providing additional information for validation and identification of patients who require further examination.



INNOVATIVE DESIGN SIMPLY INDISPENSABLE

Innovative design makes this real-time system an ideal companion for any point-of-care veterinarian. The convenience of B-mode, M-mode, Doppler and simple to access presets, makes SONIMAGE P3 simply indispensable.

Easy Access to USB port and headphone jack, touch screen stylus and 4GB microSD data card (included) which enables you to store 10,000+ images.

SONIMAGE VIEW
Online imaging and Offline review and Storage.

High Resolution Touchscreen
for easy, quick adjustment of features and functions of the device.

Versatile Transducer
The S3-5 transducer is suitable for a variety of applications like abdomen, thoracic and bladder scanning.



Pregnancy Check Mass Measurement Bladder



Changeable Transducer

SONIMAGE P3


Personal, Portable, Point-of-Care Ultrasound

Improving the Quality-of-Care at the Point-of-Care

The SONIMAGE P3 gives veterinary providers a cost-effective tool in which to offer services that require the use of ultrasound. More can now be done in-house, saving patients time while improving the level of care.

FEATURES

Weight	<ul style="list-style-type: none">• 388 grams (13.7oz)
Dimensions	<ul style="list-style-type: none">• Display unit: 11.5cm H x 6.8cm W x 2.0cm D• Display Screen: 6.85cm Hi-res colour LCD• Probe: 15cm L x ~ 6cm W x ~3.2cm D• Connecting cable: 80.0cm L
Transducer	<ul style="list-style-type: none">• S3-5 Sector Annular (3-5 MHz) – B/M Mode + PW Doppler
User interface	<ul style="list-style-type: none">• Touch screen and scroll wheel• Onscreen alphanumeric keyboard• Text and Voice Recording• Optional tablet display• PictBridge compatible printing
Measurement tools	<ul style="list-style-type: none">• Auto Volume Measurements• Linear distance calipers• Ellipse & polygon area/circumference• Beats-per-minute (BPM)
Image storage	<ul style="list-style-type: none">• 4 GB microSD data card (included)• Urgent Mode internal memory• USB connection for downloads
Software	<ul style="list-style-type: none">• SONIMAGE VIEW PC/tablet software• Backup images/data• Search and manage database• Automatic software updates
Battery	<ul style="list-style-type: none">• 75% charge in 2 hours• 100% charge in 3 hours• Lasts up to 12 hours (based on four 1-minute scans/hour)• Over 120 minutes of continuous scanning

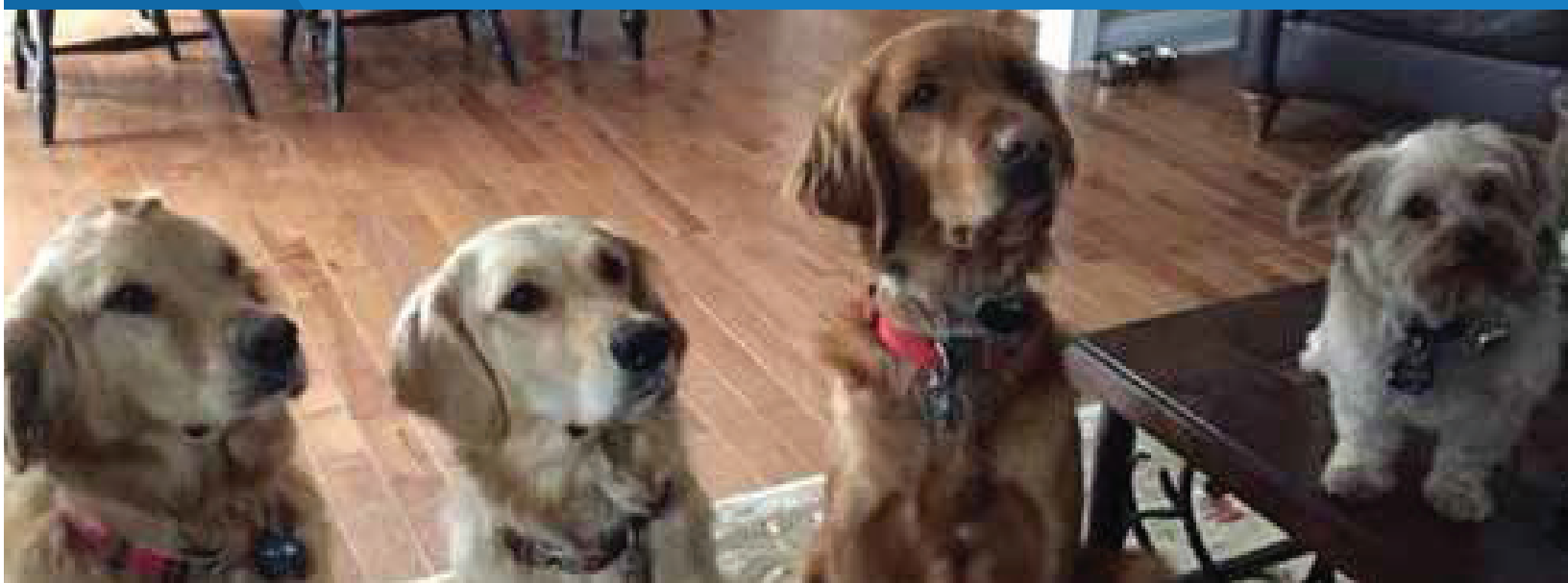

KONICA MINOLTA

Konica Minolta Medical Imaging

411 Newark Pompton Turnpike
Wayne, NJ 07470

Tel: (973) 633-1500 Fax: (973) 523-7408
konicaminolta.com/medicalusa